
# HARMONIQ

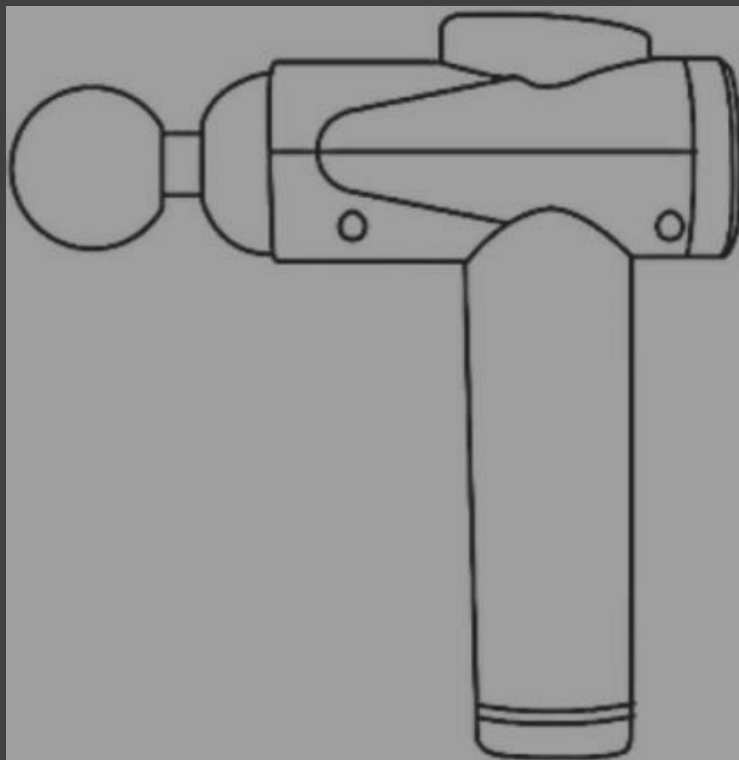User Manual | Manual de Utilizare | Használati utasítás |  
Ръководство за употреба

**English**

03-05

-

**Română**

06-08

-

**Magyar**

09-11

-

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12-14


## 1.0 OVERVIEW

Congratulations on your selection of the HARMONIQ Professional Massage Gun. HARMONIQ is family-built electronics & wellbeing brand, introducing the ultimate wellness experience with our brand's top-of-line massage gun. HARMONIQ massage guns are designed for the ones who value quality and innovation. The deep muscle fascia massage gun is a vibration massage gun which can effectively reduce the muscle pain caused by the accumulation of lactic acid from exercise and relieve the muscle discomfort of people who have been sitting for a long time at the computer. Robust materials, lavish finishes, bold geometry and outstanding performance are celebrated attributes of the HARMONIQ brand.

## 2.0 WHAT'S IN THE BOX

- Massage Gun
- 8 Interchangeable Massage Heads
- Charging Cable
- Case
- Warranty & Customer Support Card

## 3.0 PRODUCT DESCRIPTION

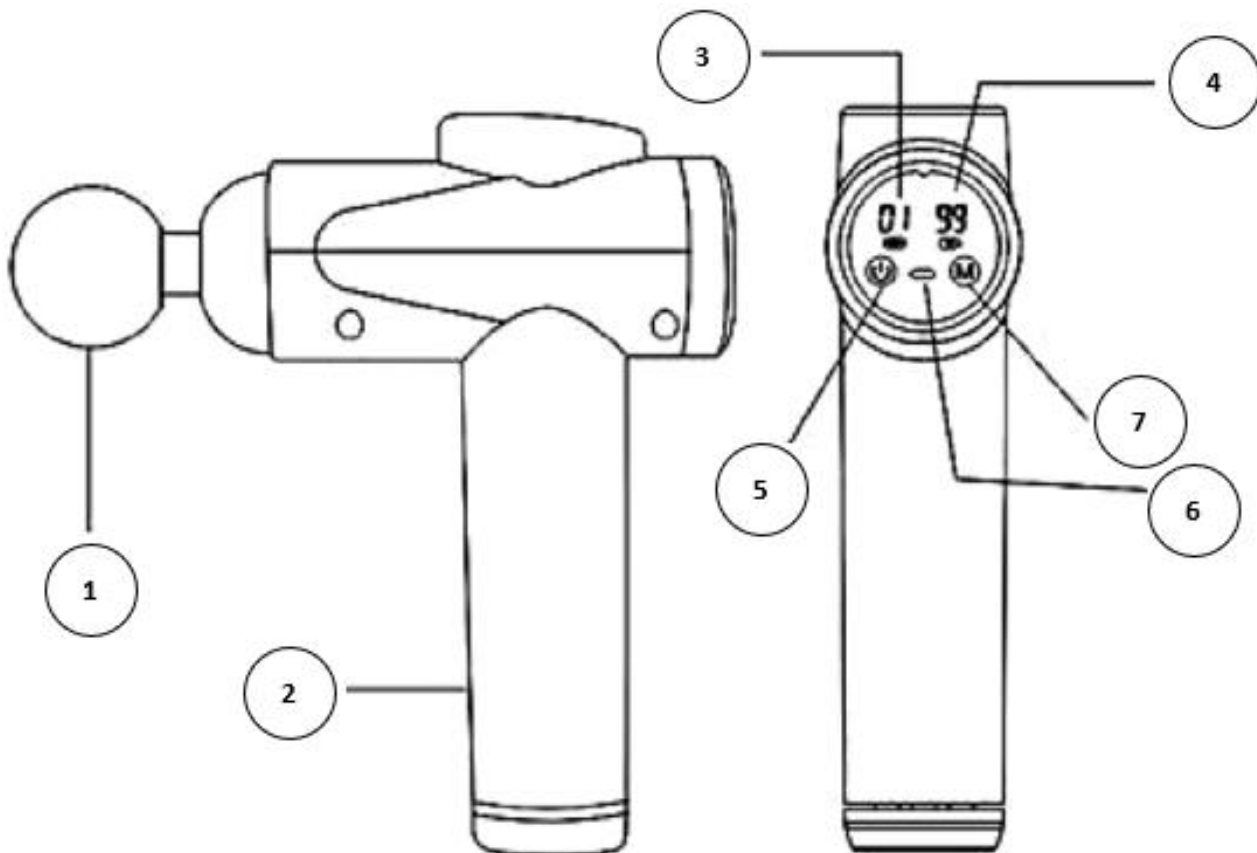1. Detachable massage head
2. Ergonomic grid
3. Strength level display
4. Battery display
5. Power button
6. Charging port
7. Mode button

#### 4.0. HOW TO USE

**Turn on:** Press and hold the Power button for 2 seconds and the massage gun will start working. The Mode button is for strength level adjustment from level 01 to 99. The first two digits on the display indicate the strength level and the last two digits indicate the battery level.

**Strength level shift:** After turning on the massage gun, the default strength level is 01. Press and hold the Mode button to quickly increase and change the level of strength until level 99.

**Turn off:** When you wish to turn off the massage gun, press and hold the power button.

The massage gun automatically protects and shuts down after 10 minutes of operation. If you wish to use it again, please press the power button again to turn it on.

#### 5.0. BATTERY INSTRUCTIONS

The last two digits on the display show the current battery level. When fully charged, it is 99. When in use, the display changes based on battery consumption, with a minimum battery change of 20. When 00 is displayed, the massage gun is low on battery. Please charge it in a timely manner.

#### WARNING:

- Please use our specially equipped USB cable with 5V power adapter to charge the product.
- Do not use a charger (adapter) with a damaged casing or an incorrect voltage to charge the product.
- Do not touch the product or charger (adapter) with wet hands to avoid product short-circuit, malfunction or electric shock. Please stop using the adapter when it is exposed to rain, liquid or severely damp.
- Do not expose the battery to high temperature or around the heating equipment, such as sunshine, heat, etc.
- Do not disassemble or refit batteries, insert foreign matters or immerse them in water or other liquids, so as to avoid leakage, overheat, fire or explosion.
- Before abandoning appliances, batteries must be removed from appliances and disposed of safely.

## 6.0 INDICATIONS

- Pain and spasm due to muscle contusion, sprain and strain
- Helps with edema fluid flow in the swollen area
- Loose and thickened connective tissue and fascia
- Reduce the accumulation of lactic acid in the muscles
- Eliminate muscle fatigue

## 7.0 CONTRAINDICATIONS

- Aneurysms, bleeding, use of blood thinners
- Heart disease, pacemaker or defibrillator
- Strain back disease
- Pregnancy, cancer
- Skin trauma, eczema
- Fracture, osteoporosis, joint disease
- Sensitive areas: head, face, cervical spine
- Tibia, dorsum of foot, dorsum of hand

In case you feel any discomfort, stop using the massage gun and consult your doctor.

## 8.0. PRODUCT PARAMETERS

Rated voltage	7.4V
Charging current	2A
Rated power	18W
Charging voltage	5V
Executive standard	GB4706.1-2005 GB4706.10-2008

In case you have questions do not hesitate to contact HARMONIQ Customer Support via scanning the QR-Code on the Warranty Card.

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